


EDUCATIONAL SCIENCE

LABORATORY TABLES

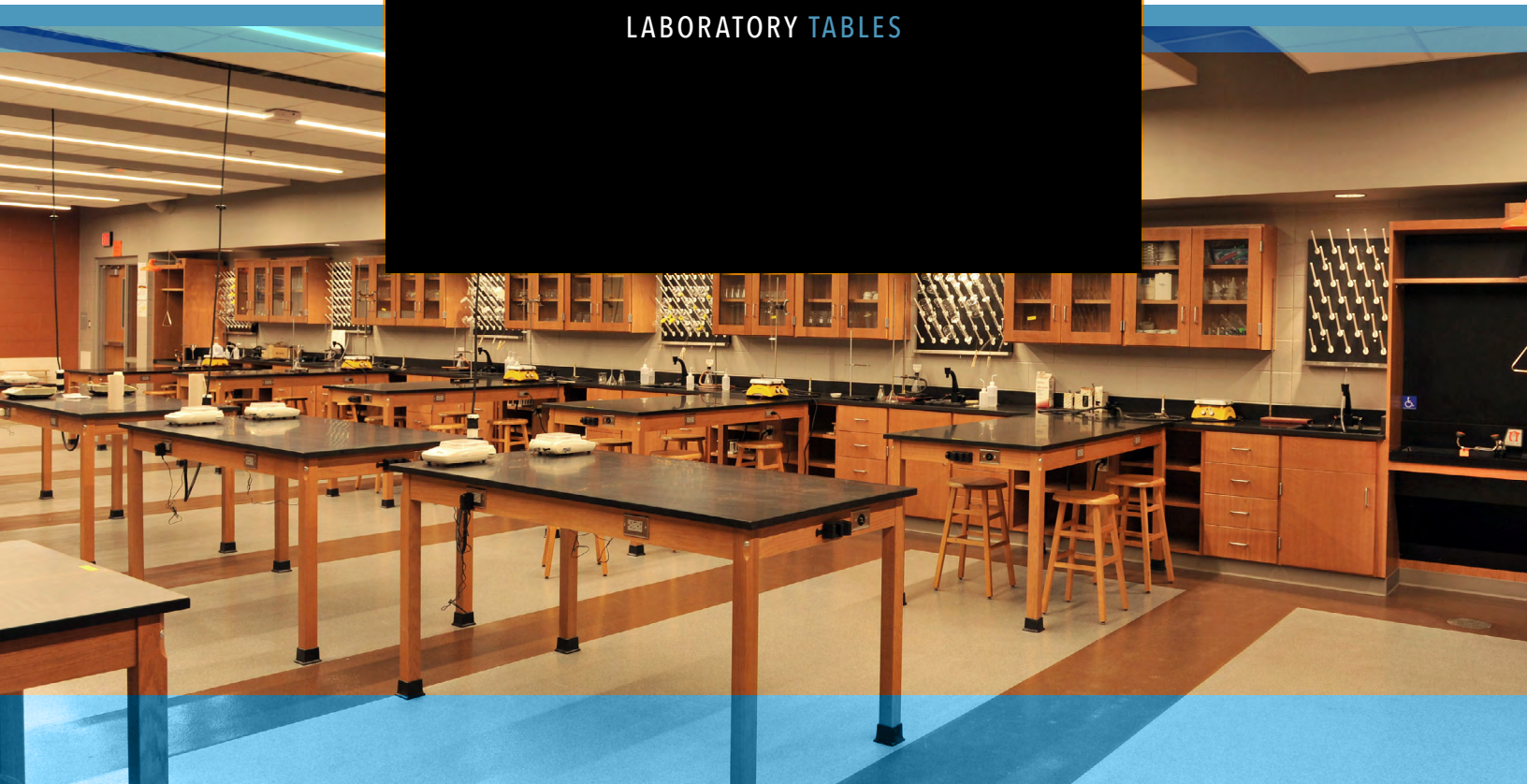

EDUCATIONAL

LABORATORY TABLES

CONFIGURATIONS

P7601K30SC

Apron Style Codes	P = PLAIN	C = BOOK COMPARTMENT	D = DRAWER	
Table Size Codes	710 = 24" X 48" 711 = 42" X 48" 712 = 30" X 48" 713 = 30" X 54" 714 = 30" X 60" 715 = 30" X 72"	716 = 21" X 48" 717 = 36" X 48" 718 = 21" X 54" 719 = 36" X 54" 720 = 24" X 54" 721 = 21" X 60"	722 = 36" X 60" 723 = 21" X 72" 724 = 36" X 72" 730 = 24" X 72" 740 = 42" X 72" 748 = 48" ROUND	760 = 24" X 60" 780 = 42" X 54" 790 = 42" X 60"
Table Surface Codes	1 = 1-1/4" PLASTIC LAMINATE 2 = 1-1/4" CHEMGUARD 4 = 1" PHENOLIC RESIN 5 = 1-3/4" MAPLE 6 = 1" EPOXY RESIN			
Wood Specie Codes	K = OAK M = MAPLE			
Table Height Codes	30 = 30" H 34 = 34" H 36 = 36" H			
Table Option Codes	N = NONE L = LEG BRACE (REQUIRED FOR ALL 42"D TABLES) E = ELECTRICAL CUTOUT AND PLATE ONLY S = LEG BRACE & STRETCHER		LE = LEG BRACE & ELECTRICAL SE = LEG BRACE, STRETCHER, ELECTRICAL SC = LEG BRACE, STRETCHER, CASTER	

Required for Tables Over 72"W // And for Some Epoxy Topped Tables.

Plain Apron Table

Are available in a variety of sizes, with several different top materials. These tables are built for years of dependable service. Legs and aprons are constructed of solid red oak, and are smoothly machined. They are finished with a chemical resistant, earth-friendly UV finish. Aprons are constructed of solid red oak and are 13/16" D x 4-1/2" H with heavy-duty 13-gauge steel corner braces. 2-1/4" square solid oak legs, with 5/16" hanger bolts installed, are attached to apron with nuts and washers supplied. All tables are equipped with non-skid adjustable glides for precise leveling and for prevention of rocking on uneven floors. Tables are supplied with black rubber leg boots which offer protection against scuffing and chemical spills. Standard table height is 30" with a weight capacity of 500 lbs.



SIZE	P. LAM	CHEMGUARD	SOLID PHENOLIC	SOLID EPOXY	NO TOP
21" x 48"	P7161K30N	P7162K30N	P7164K30N	P7166K30N	P7160K30N
21" x 54"	P7181K30N	P7182K30N	P7184K30N	P7186K30N	P7180K30N
21" x 60"	P7211K30N	P7212K30N	P7214K30N	P7216K30N	P7210K30N
21" x 72"	P7231K30N	P7232K30N	P7234K30N	P7236K30N	P7230K30N
24" x 48"	P7101K30N*	P7102K30N*	P7104K30N	P7106K30N*	P7100K30N
24" x 54"	P7201K30N*	P7202K30N*	P7204K30N	P7206K30N*	P7200K30N
24" x 60"	P7601K30N*	P7602K30N*	P7604K30N	P7606K30N*	P7600K30N
24" x 72"	P7301K30N*	P7302K30N*	P7304K30N	P7306K30N*	P7300K30N
30" x 48"	P7121K30N	P7122K30N	P7124K30N	P7126K30N	P7120K30N
30" x 54"	P7131K30N	P7132K30N	P7134K30N	P7136K30N	P7130K30N
30" x 60"	P7141K30N*	P7142K30N*	P7144K30N	P7146K30N*	P7140K30N
30" x 72"	P7151K30N	P7152K30N*	P7154K30N	P7156K30N	P7150K30N
36" x 48"	P7171K30N	P7172K30N	P7174K30N	P7176K30L	P7170K30N
36" x 54"	P7191K30N	P7192K30N	P7194K30N	P7196K30L	P7190K30N
36" x 60"	P7221K30N	P7222K30N	P7224K30N	P7226K30L	P7220K30N
36" x 72"	P7241K30N	P7242K30N	P7244K30N	P7246K30L	P7240K30N
42" x 48"	P7111K30L	P7112K30L	P7114K30L	P7116K30L	P7110K30L
42" x 54"	P7801K30L	P7802K30L	P7804K30L	P7806K30L	P7800K30L
42" x 60"	P7901K30L	P7902K30L	P7904K30L	P7906K30L	P7900K30L
42" x 72"	P7401K30L	P7402K30L	P7404K30L	P7406K30L	P7400K30L
48" ROUND	P7481K30N	P7482K30N	P7484K30N	P7486K30N	P7480K30N




SHOWN **Book Compartment Tables**

> Trapezoid Table

This unique shaped science table allows for flexibility. Use as a two student table or have two or more tables together for group work. The sturdy construction ensures that moving the tables won't cause problems. The tables are made of solid oak or maple aprons, 2-1/4" solid legs and have a steel frame supporting the apron. The wood is finished with a chemical resistant, non-emitting UV finish. The tops are your choice of 1-1/4" plastic laminate or ChemGuard, with a height of 30" with levelers or 36"H with casters.

ITEM #	WOOD	TOP	CASTERS	SIZE
TZ7141K30N	Oak	Plastic Laminate	No	30" x 60" x 30"
TZ7141M30N	Maple	Plastic Laminate	No	30" x 60" x 30"
TZ7142K30N	Oak	ChemGuard	No	30" x 60" x 30"
TZ7142M30N	Maple	ChemGuard	No	30" x 60" x 30"
TZ7601K30N	Oak	Plastic Laminate	No	24" x 60" x 30"
TZ7601M30N	Maple	Plastic Laminate	No	24" x 60" x 30"
TZ7602K30N	Oak	ChemGuard	No	24" x 60" x 30"
TZ7602M30N	Maple	ChemGuard	No	24" x 60" x 30"
TZ7141K36C	Oak	Plastic Laminate	Yes	30" x 60" x 36"
TZ7141M36C	Maple	Plastic Laminate	Yes	30" x 60" x 36"
TZ7142K36C	Oak	Chemguard	Yes	30" x 60" x 36"
TZ7142M36C	Maple	Chemguard	Yes	30" x 60" x 36"
TZ7601K36C	Oak	Plastic Laminate	Yes	24" x 60" x 36"
TZ7601M36C	Maple	Plastic Laminate	Yes	24" x 60" x 36"
TZ7602K36C	Oak	Chemguard	Yes	24" x 60" x 36"
TZ7602M36C	Maple	Chemguard	Yes	24" x 60" x 36"

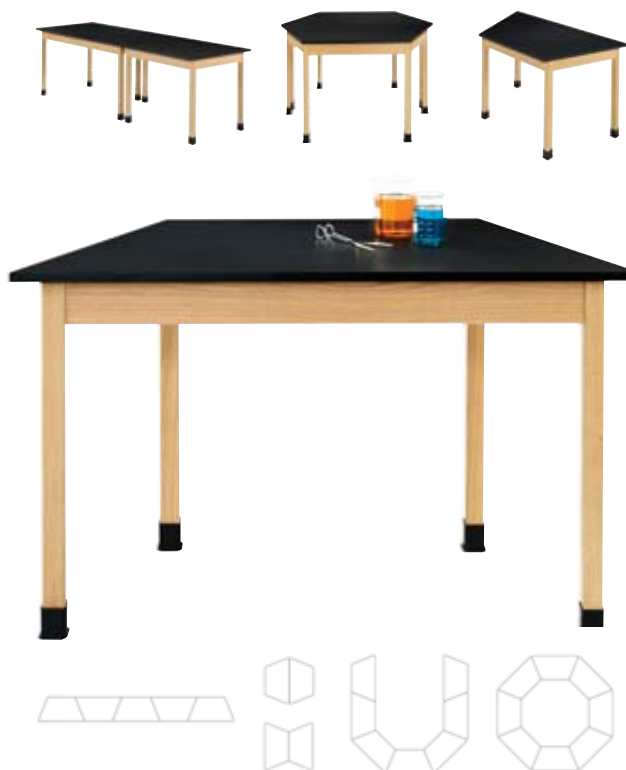

> Fixed Height Riser Table

Perfect table for STEM projects. We've taken our popular oak and maple science tables and enhanced them by adding a riser. The table features solid oak or maple (your choice) aprons and riser supports finished with a chemical resistant finish. The riser provides 10" of clearance above the work surface, is 12"D and made of the same material as the top you choose. Top options include: plastic laminate, chemical resistant laminate, phenolic resin and epoxy resin. All 72"W tables include leg braces and a stretcher for added strength. Tables stand 30"H. The top of the riser is 41" from the floor. The riser top is fully assembled, with the legs needed to be attached.

SIZE	PLASTIC LAMINATE OAK	PLASTIC LAMINATE MAPLE	CHEMGUARD OAK	CHEMGUARD MAPLE	PHENOLIC OAK	PHENOLIC MAPLE	EPOXY OAK	EPOXY MAPLE
30" x 54"	PR7131K30N	PR7131M30N	PR7132K30N	PR7132M30N	PR7134K30N	PR7134M30N	PR7136K30N	PR7136M30N
30" x 60"	PR7141K30N	PR7141M30N	PR7142K30N	PR7142M30N	PR7144K30N	PR7144M30N	PR7146K30N	PR7146M30N
30" x 72"	PR7151K30S	PR7151M30S	PR7152K30S	PR7152M30S	PR7154K30S	PR7154M30S	PR7156K30S	PR7156M30S



> Adjustable Wooden Leg Table

The A7000 table has the same construction as our standard table with the exception of adjustable legs. The aprons are solid oak and are available in four types of tops: epoxy, phenolic resin, plastic laminate, and ChemGuard. Legs can adjust from sitting height of 30" to standing height of 36" by removing the bolt on the leg surface. The table is finished with a chemical resistant finish.

ITEM #	TOP	SIZE
A7101K	Plastic Laminate	48" x 24" x 30"-36"
A7102K	ChemGuard	48" x 24" x 30"-36"
A7104K	Phenolic	48" x 24" x 30"-36"
A7106K	Epoxy	48" x 24" x 30"-36"
A7141K	Plastic Laminate	60" x 30" x 30"-36"
A7142K	ChemGuard	60" x 30" x 30"-36"
A7144K	Phenolic	60" x 30" x 30"-36"
A7146K	Epoxy	60" x 30" x 30"-36"
A7201K	Plastic Laminate	54" x 24" x 30"-36"
A7202K	ChemGuard	54" x 24" x 30"-36"
A7204K	Phenolic	54" x 24" x 30"-36"
A7206K	Epoxy	54" x 24" x 30"-36"
A7301K	Plastic Laminate	72" x 24" x 30"-36"
A7302K	ChemGuard	72" x 24" x 30"-36"
A7304K	Phenolic	72" x 24" x 30"-36"
A7306K	Epoxy	72" x 24" x 30"-36"
A7601K	Plastic Laminate	60" x 24" x 30"-36"
A7602K	ChemGuard	60" x 24" x 30"-36"
A7604K	Phenolic	60" x 24" x 30"-36"
A7606K	Epoxy	60" x 24" x 30"-36"



With Stretchers, Leg Braces & Casters

> Book Compartment Tables

Are built with the same quality techniques as the plain apron table. They are finished with a chemical resistant, earth-friendly UV finish. These tables have either two or four book compartments (depending on size) installed in the apron. The book compartments measure 14" W x 13-5/8"D x 3-7/8"H. The table has a weight capacity of 500 lbs. On tables 72"W or over, we recommend ordering with leg braces and stretchers. (Shipped KD with tops attached)

SIZE	COMP.	PLASTIC LAMINATE	CHEMGUARD	SOLID PHENOLIC	SOLID EPOXY	NO TOP
21" x 48"	2	C7161K30N	C7162K30N	C7164K30N	C7166K30N	C7160K30N
21" x 54"	2	C7181K30N	C7182K30N	C7184K30N	C7186K30N	C7180K30N
21" x 60"	2	C7211K30N	C7212K30N	C7214K30N	C7216K30N	C7210K30N
21" x 72"	2	C7231K30N	C7232K30N	C7234K30N	C7236K30N	C7230K30N
24" x 48"	2	C7101K30N*	C7102K30N*	C7104K30N	C7106K30N*	C7100K30N
24" x 54"	2	C7201K30N*	C7202K30N*	C7204K30N	C7206K30N*	C7200K30N
24" x 60"	2	C7601K30N*	C7602K30N*	C7604K30N	C7606K30N*	C7600K30N
24" x 72"	2	C7301K30N*	C7302K30N*	C7304K30N	C7306K30N*	C7300K30N
30" x 48"	2	C7121K30N	C7122K30N	C7124K30N	C7126K30N	C7120K30N
30" x 54"	2	C7131K30N	C7132K30N	C7134K30N	C7136K30N	C7130K30N
30" x 60"	2	C7141K30N*	C7142K30N*	C7144K30N	C7146K30N*	C7140K30N
30" x 72"	2	C7151K30N	C7152K30N*	C7154K30N	C7156K30N	C7150K30N
36" x 48"	2	C7171K30N	C7172K30N	C7174K30N	C7176K30L	C7170K30N
36" x 54"	2	C7191K30N	C7192K30N	C7194K30N	C7196K30L	C7190K30N
36" x 60"	2	C7221K30N	C7222K30N	C7224K30N	C7226K30L	C7220K30N
36" x 72"	2	C7241K30N	C7242K30N	C7244K30N	C7246K30L	C7240K30N
42" x 48"	4	C7111K30L	C7112K30L	C7114K30L	C7116K30L	C7110K30L
42" x 54"	4	C7801K30L	C7802K30L	C7804K30L	C7806K30L	C7800K30L
42" x 60"	4	C7901K30L	C7902K30L	C7904K30L	C7906K30L	C7900K30L
42" x 72"	4	C7401K30L	C7402K30L	C7404K30L	C7406K30L	C7400K30L

TABLE PARTS

UNIT #	DESCRIPTION	UNIT #	DESCRIPTION
L30-REP*	30"H Leg (1 each)	L3034	34"H Leg (set of 4)
L34-REP	34"H Leg (1 each)	L3036*	36"H Leg (set of 4)
L36-REP*	36"H Leg (1 each)	100076X4*	Leg Boots (set of 4)
L3030*	30"H Leg (set of 4)	100728X4*	Floor Mounting Brackets (set of 4)



> Adjustable Height Riser Table

This table has everything you could want. It adjusts easily from 27-3/4" to 39" (42"H with 3" spacer provided) in 2 cm increments in less than 30 seconds without the use of tools. Apron is solid oak and finished with a chemical resistant finish. It has a riser that provides extra workspace for each to put tools or laptop. The table has three top options: plastic laminate, chemical resistant laminate and phenolic resin. The riser's horizontal surface matches your choice of top surface. The 12"D riser is supported by six 1" diameter powder coated steel tubes. The riser is 10" above the table surface. Tabletop measures 30"D x 60"W.

ITEM #	TOP
PR8141K	Plastic Laminate
PR8142K	Chemical Resistant Laminate
PR8144K	Phenolic Resin



> Adjustable Height Table

A perfect solution for ADA requirements that can adjust in 2 cm increments from 27-3/4"H to 39" (42"H with 3" spacer provided). Can also be adjusted for use as an incline surface for physics. Solid 4-1/2" oak aprons are attached to your choice of three top options: plastic laminate, ChemGuard, and phenolic. Note: phenolic resin tops are heavier and may require two people to adjust the height. The wood is finished with a chemical resistant, earth-friendly UV finish. The table has a weight capacity of 500 lbs. Optional casters sold separately.



Shown with Optional Castors

SIZE	PLASTIC LAMINATE	CHEMGUARD	SOLID PHENOLIC
24" x 48"	P8101K*	P8102K*	P8104K
24" x 54"	P8201K*	P8202K*	P8204K
24" x 60"	P8601K*	P8602K*	P8604K*
24" x 72"	P8301K*	P8302K*	P8304K
30" x 60"	P8141K*	P8142K*	P8144K
Castors	100709X4*	(Set of 4)	

> **Storage Cabinet Tables**

Are constructed of solid oak and hardwood veneers, this table features a center cabinet with lock and rubber base molding. This table is available with plain or book compartment aprons and comes in both single (2-student) and double-faced (4-student) models. The center cabinet is available with two locking drawers or as a plain base cabinet with adjustable shelf. These tables are designed for use at sitting height. They are finished with a chemical resistant, earth-friendly UV finish. This unit is available with our complete line of top materials to satisfy every requirement or budget.



2-STUDENT SINGLE	WEIGHT	4-STUDENT DOUBLE	WEIGHT	APRON	CABINET	TOP
2101K	265 lbs.	2301K	380 lbs.	Plain	Plain	Plastic Laminate
C2101K	275 lbs.	C2301K	400 lbs.	Compartment	Plain	Plastic Laminate
2201K	275 lbs.	2401K	400 lbs.	Plain	Drawers	Plastic Laminate
C2201K	275 lbs.	C2401K	410 lbs.	Compartment	Drawers	Plastic Laminate
2102K	265 lbs.	2302K	380 lbs.	Plain	Plain	ChemGuard
C2102K	275 lbs.	C2302K	400 lbs.	Compartment	Plain	ChemGuard
2202K	275 lbs.	2402K	400 lbs.	Plain	Drawers	ChemGuard
C2202K	275 lbs.	C2402K	410 lbs.	Compartment	Drawers	ChemGuard
2104K	316 lbs.	2304K	443 lbs.	Plain	Plain	Solid Phenolic
C2104K	326 lbs.	C2304K	463 lbs.	Compartment	Plain	Solid Phenolic
2204K	330 lbs.	2404K	463 lbs.	Plain	Drawers	Solid Phenolic
C2204K	335 lbs.	C2404K	468 lbs.	Compartment	Drawers	Solid Phenolic
2106K	340 lbs.	2306K	500 lbs.	Plain	Plain	Solid Epoxy
C2106K	350 lbs.	C2306K	520 lbs.	Compartment	Plain	Solid Epoxy
2206K	350 lbs.	2406K	520 lbs.	Plain	Drawers	Solid Epoxy
C2206K	355 lbs.	C2406K	525 lbs.	Compartment	Drawers	Solid Epoxy



SINGLE FACE UNIT (66"W X 24"D X 30"H)

- 1 Single-Faced Cabinet (24"W x 22"D x 29"H)
- 2 Book Compartments (12"W x 13"D x 3"H)
- 2200 Model: 2 Drawers Unit (11-1/2"W x 18"D x 5-1/2"H)
- 2100 Model: Cabinet Base Unit w/full door have an adjustable shelf, with a 40 lb. capacity per sq. foot.



DOUBLE FACE UNIT (66"W X 42"D X 30"H)

- 1 Double-Faced Cabinet (24"W x 40"D x 29"H)
- 4 Book Compartments Unit (12"W x 13"D x 3"H)
- 2400 Model: 4 Drawer Unit (11-1/2"W x 18"D x 5-1/2"H)
- 2300 Model: Cabinet Base Unit w/full door have an adjustable shelf, with a 40 lb. capacity per sq. foot.



Back of Cabinet

> Add-A-Cabinet

This unique table gives you the ability to change it as your needs change. Basic unit comes with a 60"W x 30"D x 36"H table that is reinforced with welded steel angle understructure. Cabinets can be mounted to the underside of the table and attached using a solid bolt connection. The unit comes standard with one cabinet, but there is enough room for two cabinets to be hung from the table. The cabinets measure 24"W x 22-1/2"D x 28-1/2"H. Choose from three cabinet options: drawer over door (with one adjustable shelf), Single door (with one adjustable shelf) and Drawer (two upper drawers and a file drawer). The doors and drawers lock. The top is 1" solid epoxy resin. The unit is finished with a non-emitting UV finish.

ITEM #	PRODUCT
H7146K36S-19	Table w/Drawer over Door Cabinet
H7146K36S-25	Table w/Drawer Cabinet
H7146K36S-33	Table w/Single Door Cabinet



M19-2422-H30K (Shown w/30"H Table)

> "M" Series Cabinets

Perfect to use as a standalone unit, or to fit underneath our best-selling tables. They are finished with a chemical resistant, earth-friendly UV finish. The cabinets roll on 4" heavy-duty ball-bearing casters with brakes. Available for sitting (24"H units) or standing (30"H units) tables. Shipped assembled except for the casters. (Table sold separately)



M96-2422-H24K



M40-2422-H30K



M33-2422-H30K

ITEM #	FOR TABLE HGT.	LAYOUT	TOP	SIZE
M95-2422-H24K	30"	1 Door/Hinged Right	Veneered	22"D x 24"W x 24"H
M96-2422-H24K	30"	1 Door/Hinged Left	Veneered	22"D x 24"W x 24"H
M95-2422-H30K	36"	1 Door/Hinged Right	Veneered	22"D x 24"W x 30"H
M96-2422-H30K	36"	1 Door/Hinged Left	Veneered	22"D x 24"W x 30"H
M18-2422-H30K	36"	1 Door/Drawer, Hinged Right	Veneered	22"D x 24"W x 30"H
M19-2422-H30K	36"	1 Door/Drawer, Hinged Left	Veneered	22"D x 24"W x 30"H
M40-2422-H30K	36"	2 Drawer	Veneered	22"D x 24"W x 30"H
M33-2422-H30K	36"	3 Drawer	Veneered	22"D x 24"W x 30"H



3216K (Shown w/Two 36"H Plain Apron Tables)



3216K



3226K



3236K



3246K

> Wall Service Bench

Need all your services in one place or looking to add a wall of cabinets that include services? The Wall Service Bench offers storage and services all in one unit. Constructed of solid hardwood and oak veneer plywood. The Wall Service Bench features your choice of a 1" solid epoxy resin or phenolic resin top with a 18"W x 15"D x 11"H epoxy resin sink and two combination cold water/gas fixtures. The sink area allows access to the plumbing and provides a large storage space with two GFI protected AC duplex electrical outlets. There is enough room in the sink area to allow two 24"D tables to be butted against the unit for use as a peninsula work area. The outer ends of the unit offer four storage solutions:

- Doors with Adjustable Shelves
- Doors/Drawer Combination w/1 Adjustable Shelf
- 4 Drawers/Doors w/1 Adjustable Shelf
- 4 Drawers

All drawers and doors have locks and ADA compliant pulls. The unit is finished with a chemical resistant, earth-friendly UV finish and includes rubber base molding. Dimensions: 108"L x 24"D x 36"H.

EPOXY	PHENOLIC	CABINET STYLE
3216K	3214K	Door
3226K	3224K	Door/Drawer
3236K	3234K	Door/4 Drawer
3246K	3244K	Drawer

Note: Units are built to be installed against a wall and do not have finished backs.